
Aprilla

EN

FR

TR

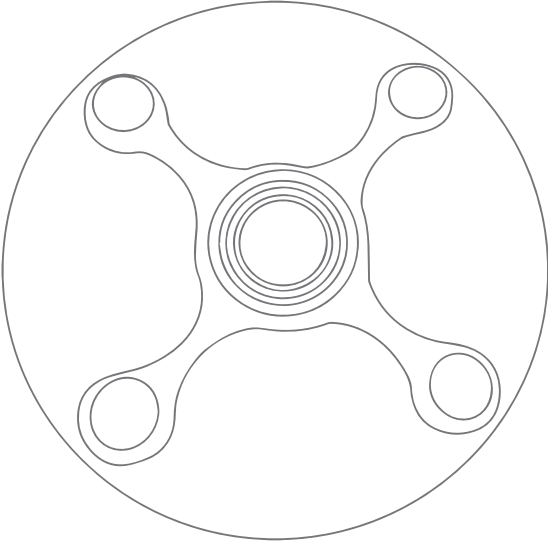ABS 1030

Dijital Banyo Baskülü

Digital Bathroom Scale

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user manual

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- *This device is designed and manufactured for home use. It is not suitable for commercial or industrial purposes.*
- *This device may be interest of children, but as you know it is not a toy. In this reason, device being where children can't reach place.*
- *Other information in this guide to see the section "**Safety Instructions**".*
- *Please read the manual, you can start using the product after that.*

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Children Being Supervised not to Play with the Appliance

this appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.

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Hello,

First of all thank you for choosing Kiwi product. Please read this manual carefully before using your device. Please keep the manual as it may need in future. This guide is designed to help you to achieve maximum performance.

We wish you a pleasant use.

Aprilla

SPECIFICATIONS

- Maximum weight capacity: 150kg/330lb/23.6st
- Minimum measurable weight: 2.5kg/5.5lb/0.4st
- Graduation: 0.1KG/0.21lb
- User selectable units: kg/lb/st
- Low-battery indicator: LO
- Overload indicator: EEEE
- Operating temperature: 5°C ~ 40°C

EN

INSTRUCTIONS FOR USE

1. Get one CR2032 Lithium battery put in the battery seat of the back of scale
2. Choice the weight unit and move to the right position.
3. The scale must put on the smooth hard floor.
4. Turn on scale by standing on the scale directly
5. The display will show 0.0kg
6. Always stand in the center of the platform, thus enabling the scale to take your weight exactly.
7. Within seconds, personal weight is established electronically and displayed
8. Once you step off the scale after weighing, it will switch itself off automatically after few seconds.

TROUBLE SHOOTING

- Display "EEEE" .if weight exceeds the maximum capacity of 150kg/330lb/23.6st. "EEEE" Appears on the display, Please step off the scale immediately, otherwise the scale may be damaged.
- Display "Lo". Low battery. The user is advised to replace a pair of new battery.

PRECAUTIONS

- Avoid overload, other wise, permanent damage to the scale will occur.
- Do not disassemble the product other than replacing the battery.
- Clean after use with a lightly damper cloth, do not use solvents or immerse the products in water.
- Avoid excessive impact or vibration to the unit, such as dropping the product onto the floor.

NOTICE

- Keep the glass platform is dry for fear body will be slip. Note the scale guard against wet.
- Do not put the scale on the soft floor (like on the carpet), otherwise the measure is inaccurate.
- Do not put any other things on the glass platform if you not use.

IMPORTANT: Remove the plastic under the battery.

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Correct Disposal of this product

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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© 2016 BMVA. All rights reserved. Aprilla logo is owned by BMVA and registered. All other trademarks and signs are property of their owners. BMVA is not responsible for printing and composing mistakes.