
Aprilla

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Body Fat Analysis Digital Bathroom Scale

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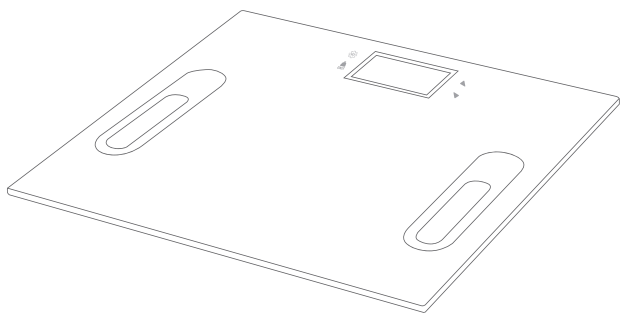- *This device is designed and manufactured for home use. It is not suitable for commercial or industrial purposes.*
- *This device may be interest of children, but as you know it is not a toy. In this reason, device being where children can't reach place.*
- *Other information in this guide to see the section "**Safety Instructions**".*
- *Please read the manual, you can start using the product after that.*

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Children Being Supervised not to Play with the Appliance

This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliances by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.

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Hello,

First of all thank you for choosing Aprilla product. Please read this manual carefully before using your device. Please keep the manual as it may need in future. This guide is designed to help you to achieve maximum performance.

We wish you a pleasant use.

Aprilla

GENERAL SAFETY INSTRUCTIONS

Read the operating instructions carefully before putting the appliance into operation and keep the instructions including the warranty, the receipt and, if possible the box with the in-ternal packing. If you give this device to other people, please also pass on the operating instruction.

EN

- The device is designed exclusively for private and not for commercial use.
- Only use the device for the intended purpose.
- Do not use a damaged device.
- Keep away from heat, direct sunlight, damp (never immerse in liquids) and sharp edges.
- Do not try to repair the appliance on your own. Always contact an authorized technician.
- Pay careful attention to the following "Special Safety Instructions".

SYMBOLS IN THESE INSTRUCTIONS FOR USE

Important information for your safety is specially marked. It is essential to comply with these instructions in order to avoid accidents and prevent damage to the machine

WARNING:

This warns you of dangers to your health and indicates possible injury risks.

CAUTION

This refers to possible hazards to the machine or other objects.

NOTE: This highlights tips and information.

Children

- In order to ensure your children's safety, please keep all packaging (plastic bags, boxes, polystyrene etc.) out of their reach.

WARNING!

Caution! Do not allow small children to play with the foil as there is a danger of suffocation!

- Children should be supervised at all times in order to ensure that they do not play with device.

WARNING:

If plastic surface is moist or covered with water droplets, please wipe the surface dry as there is a danger of slipping.

CAUTION:

- Handle the balance with care, as it is a precise measuring instrument. Do not prop it or jump on it.
- Do not load the appliance with more than 180kg to avoid damaging the device.
- If the scales are carried or not going to be used, please store them only in a horizontal position.
- The body fat analysis scale is not suitable for the following persons:
 - Children under 10 years and adults over 85 years
 - Persons with symptoms of fever, oedema, or osteoporosis
 - Persons undergoing treatment for dialysis
 - Persons who use a pacemaker
 - Persons who take cardiovascular medication
 - Pregnant women
 - Sports persons who perform over 5 hours of intensive sporting activities each day
 - Persons with a heartbeat lower than 60 beats/minute

HANDLING BATTERIES**WARNING:**

- Do not expose the battery to high temperatures or direct sunlight. Never throw batteries into the fire. There is a danger of explosion!
- Keep batteries away from children. They are not a toy!
- Do not open the batteries by force.

WARNING:

- Avoid contact with metallic objects. (Rings, nails, screws et cetera.) there is a danger of short-circuiting!
- As a result of a short-circuit batteries may heat up considerably or even catch fire. This may result in burns.
- For your safety the battery poles should be covered with adhesive strips during transport.
- If a battery leaks do not rub the liquid into contact with the liquid they should be washed thoroughly and eyes rinsed out with clear water, if symptoms persist please consult a doctor.

INSERTING THE BATTERIES

Where required, the supplied batteries have been protected for transport with plastic foil.

1. Remove this foil before inserting the batteries into the appliance.
2. Open the battery compartment on the back.
3. Insert 1 pcs Lithium battery, please ensure that the polarity is correct (see the floor of the battery compartment)!
4. Close the battery compartment again.

Aprilla

5. If the device is not used for a prolonged period of time, please take the batteries out in order to Prevent leakage of battery acid.

6. If the message "LO" appears on the display, replace the batteries.


CAUTION:

- Different battery types or new and used batteries must not be used together.
- Batteries are not to be disposed of together with domestic wastes. Dispose of the used batteries at the relevant collection bin or at your dealer's.

EN

LOCATION OF CONTROLS

1. Contact surfaces
2. Display
3. DOWN-button
4. UP-button
5. SET-button
6. Unit

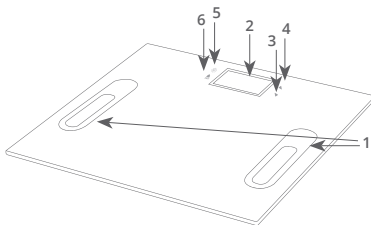Display Symbol

Kg/lb/st

Meaning

Body weight

P0-9

User memory space

Male

Female

CM(FT;In)

Body height

Age

Age

FAT

Body fat display

TBW

Body water display

%

Percentage-body fat/body water/ weight

TECHNICAL SPECIFICATIONS

- **Power supply:** Batteries 2XCR2032
- **Net weight of device:** 1.65KG
- **Memory capacity:** 10
- **Mass alignment:** from 75 to 225 cm(2'5"-7'4 1/2 "Ft:in)
- **Age:** from 10 to 85 years
- **Minimum weight:** from 2.5kg(5.5lb/0:4st)
- **Maximum weight:** 150kg(330.7lb/23:6st)
- **Weight measurement unit:** 100g increments
- **Body fat measuring unit:** 0.1%
- **Body fat measuring range:** 4.0%-50.0%
- **Body water measuring range:** 35.0%-85.0%

Subject To Technical Changes Without Prior Notice!

This device has been tested according to all relevant current CE guidelines, such as electromagnetic compatibility and low voltage directives, and has been constructed in accordance with the latest regulations.

INFORMATION ON MEASURING BODY FAT

WARNING:

Do not use the scales if you use a pacemaker or another medical implant with electronic components. The function of the implant may be impaired the measurement.

NOTE: The device can only function correctly if you stand on the contact surface with bare, dry feet! The main principle of this device is based on the measurement of the electric impedance within the human body. In order to determine the actual body data, the device sends a weak, unnoticeable electric signal through your body. This signal is generally safe! However, please read our special safety instructions! The method of measurement is referred to as "bioelectric impedance analysis" (BIA), a factor which is based on the relationship between body weight and associated body water percentages as well as on other biological data (age, gender, body weight). There are heavy persons with a well-balanced percentage between fat and muscle mass, and there are seemingly slim persons with significantly high body fat values. Fat is a factor that puts serious strain on the body, in particular the cardiovascular system. Therefore, controlling both values (fat and weight) is a prerequisite for health and fitness.

Body Fat

Body fat is an important component of the body. From protecting the joints as a vitamin repository to regulating the body temperature, fat plays a vital role in the body. Thus the goal is not radically reduce the proportion of body fat, but rather to achieve and maintain a healthy balance between muscle and fat.

Use

The device has a very sensitive sensor. If you want to adjust the settings, simply touch the buttons gently with your finger. If the device does not respond to touch, clean the glass surface with a soft, dry cloth.

NOTE: The scales are factory set with the units "kg" and "lb". To change the settings to "lb" or "kg" and "ft:in", please touch the "unit" button

EN

Weighing Function

Place the scale on a hard, flat surface (avoid carpets). An uneven floor contributes to inaccurate

Measurements.

- Stand on the scales with both feet. The device switches on automatically.
- Distribute your weight equally and wait for the scale to determine your weight.
- The digits on the display will flash before your exact weight will be shown.
- In order to extend the battery life, the device will automatically turn off 10 seconds after you have stepped off the scale.

Saving Personal Data

NOTE:

- The scale has default values that you can change, as described below.
- Carry out the following steps in sequence and within a few seconds. Otherwise, the device will automatically cancel input mode.
- Press the button SET. The memory space will flash on the display (USER0). Confirm the setting by pressing SET or select another memory space by pressing the UP/DOWN button. Then press SET to go to the next setting.
- The symbol for "male" or "female" will flash. Confirm the setting by pressing SET or use the UP/DOWN button to choose the symbol for your gender. Then press to go to the next setting.
- The height setting will flash. Select your height by pressing UP or DOWN button. Holding down one of the buttons will accelerate the process. Confirm your entry by pressing SET.
- The setting AGE will flash. Select your age by pressing UP or DOWN button. Holding down one of the buttons will accelerate the process. Confirm your entry by pressing SET.

NOTE:

- The device can only function correctly if you stand on the device with bare, dry feet!
- If necessary wipe the contact surface improve results.
- Due to the significant fluctuation of body fat, it is different to attain a standard value. on average The optimal testing time is between 6 p.m. And 8 p.m, since the body fat remains very Consistent during this period of time.
- Press the SET button. the preset stored position date appears on the display (USER flash). Use the UP/DOWN button to select the stored position at which your date is saved.
- Wait for the USER display and 0.0kg(0.0lb/st) to stop flashing.
- Stand on the contact surfaces(1) of the scales with both feet.
- Firstly your weight will flash on the display.
- After a few seconds, the following results will be shown:
 - the body fat content in%(FAT symbol)
 - the body water content in %(TBW symbol)
 - the results analysis in the form of a bar

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Correct Disposal of this product

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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This device appropriate for regulations according to Waste Electrical and Electronic Equipment (WEEE).

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© 2018 BMVA. Tüm hakları saklıdır. Aprilla logosu BMVA'ya aittir. Diğer tüm logo ve markalar sahiplerine aittir. Yazım ve baskı hatalarından BMVA sorumlu değildir. © 2018 BMVA. All rights reserved. Aprilla logo is owned by BMVA and registered. All other trademarks and signs are property of their owners. BMVA is not responsible for printing and composing mistakes.